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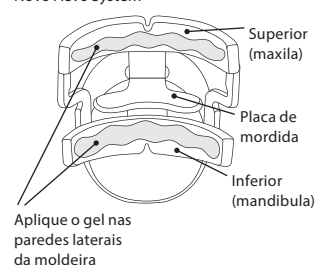
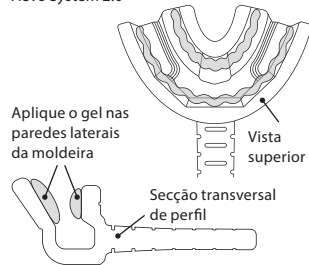
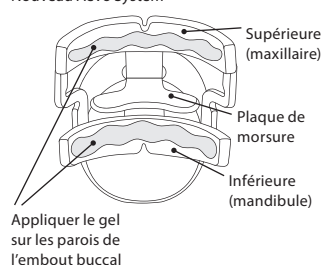
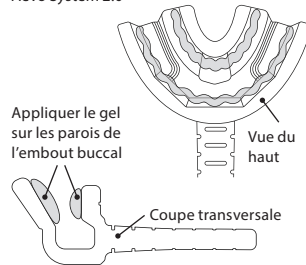
Français

SmileSonica Inc.
11514 – 168 Street NW
Edmonton, Alberta
Canada T5M 3T9

Oral Ultrasound Gel Instructions for Use

Applying Ultrasound Gel to the Mouthpiece(s)

WARNING: Ensure that the ultrasound gel is not expired before use. Only use the ultrasound gel provided with the Aevo System™. Other ultrasound coupling gels available on the market are not designed or tested for use with the Aevo System™. If any adverse reaction to the gel occurs, contact your dentist.

Before starting a treatment, tear open a corner of the gel pouch. All the gel from the pouch should be applied as a thin layer (3–4 mm thick) to the inside of the mouthpiece (see figure). This ensures contact between the mouthpiece and your gums.

NOTE: If you are using the Aevo System 2.0 Dual model, use a second gel pouch for the second mouthpiece.

Do not swallow the gel. Rinse out your mouth with water after each treatment.

NOTE: Use the ultrasound gel during every treatment or the treatment may not be effective.

IMPORTANT:

- Gel pouches are single use. This prevents contamination or infection.
- Store gel at room temperature between 15°C and 30°C.
- Expiry date format on gel pouches is YYYY-MM-DD.

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